
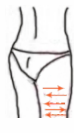


Massage Methods

Be sure to start with the lowest energy to prevent burns. increase slowly

Face (Treatment time: 30mins)	Follow the vertical direction to the ear shown by the hand gesture 	Follow the hand gesture clockwise and circle the face 	
THIGH (Treatment time: 30mins)	From down to upper pushing to the groin to dredge the lymph 	Horizontally along the thigh follow the hand gesture back and forth 	Also push by back and forth press motion to decompose fatness 
BACK (Treatment time: 20-30mins)	Lift from bottom to top 	By lymph drainage system as illustrated, pull the toxin to the lymph node 	
ARM (Treatment time: 30mins)	From the inner elbow to the armpit follow lymph drainage system as illustrated 	From the elbow to the outer, do tightening gesture as illustrated 	
ABDOMEN (Treatment time: 20-30mins)	From small to large circles, in clockwise direction, circle around the navel 	Follow the lymph drainage system, as illustrated hand gesture, transfer the toxin to around the groin region 	At the belly region pull towards the groin region as illustrated from the two sides inwards 

MINI- RF

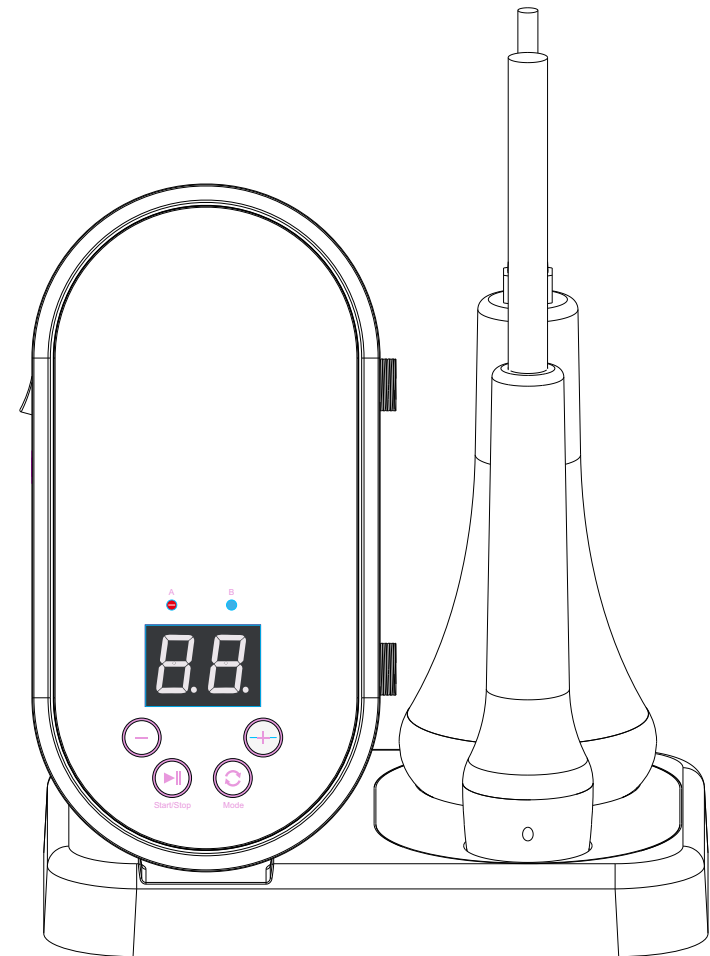

Contraindication

1. Pregnant women or women during menses.
2. Epileptic
3. Patients with malignancy.
4. Patient whose wound after operation has not healed up.
5. Acute inflammation or epidemical patients.
6. Those with heart diseases or with heart pacemaker.
7. Those with kidney (gall-stone) disease.
8. Those with embedded metal object or silica gel.
9. Those in menses, under birth control, with emiction incontinence period, or undergone belly operation.
10. Those with genetic hypersensitivity.


Product parameters

Energy output: Level 20


Power: $\leq 75W$

Product size: 13.5*18*21CM

Power input: 100--240V

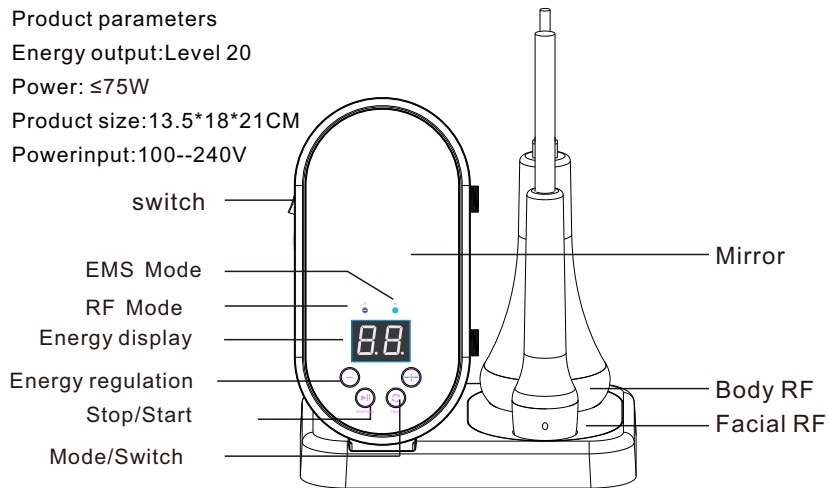

Note: When the indicator light is flashing, it indicates that the machine is working.
Red light (for RF mode). Blue light (for micro-current mode)

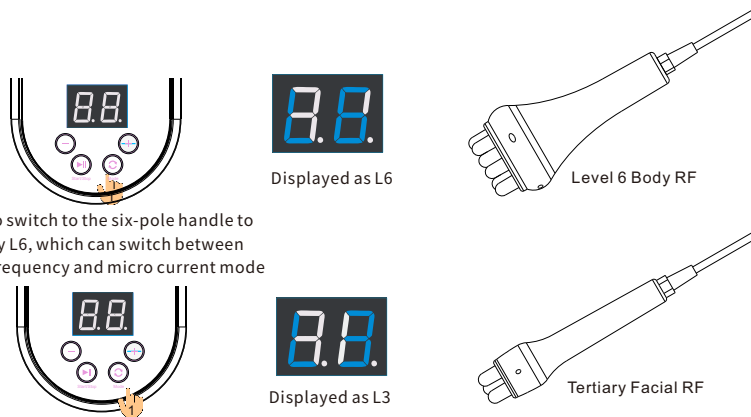

Click to switch to the six-pole handle to display L6, which can switch between radio frequency and micro current mode

Click to switch to the three-pole handle to display L3, which can switch between radio frequency and micro current mode

Function and performance

For skin tissue, When acting on the skin, the collagen fibers in the dermis layer can be heated to 45°C to 65°C, which will cause the contract of collagen fibers, leading to the proliferation of collagen so as to remove wrinkles, reduce pigmentation and restore skin elasticity.

Mode Benefits

EMS---Strengthen muscle, Reduce puffiness, Encourage lymphatic drainage, Improve skin circulation.

RADIOFREQUENCY THERAPY---Improve skin's elasticity, Tighten and firm, Improve skin's overall texture and volume.

Using this device is helpful for skin care as follows:

1. Manage wrinkles or rejuvenate the skin The device delivers high frequency radiofrequency energy to the skin, stimulating the production of collagen in the skin layers. In this case, the stimulated collagen tissue is activated, thus helping to prevent wrinkles caused by stress and helping to produce collagen in the skin layers. The stimulated collagen tissues are activated and help to speed up the rejuvenation of the skin.

2. This device relieves freckles.
This device alleviates freckles caused by deteriorating blood circulation by delivering high frequency radiofrequency energy to the skin to accelerate blood circulation. It also achieves a whitening effect by reducing the production of melanin. The antibiotic effect is activated when heat is generated in the skin tissue, which then normalises the sebaceous glands and regulates the production of sebum.

Step 1: Cleanse the skin

Step 2: Rub on the special gel

Step 3: Treat with radio frequency for 20-30 minutes

Step 4: Cleanse the skin and do a moisturizing mask

You can rub some moisturizing essence on your body.

Note: Within three days of a treatment, apply moisturizing extract to the treated area



RF energy

Skin irritation

Produces more collagen