

Warnings:Important Information For Users In order for products to perform properly,we recommend the following applicable warranties.

(1)When using the device, Do not look into the laser beam either directly or by any kind of optical equipment to prevent discomfort (2)When irradiating the treatment site, please do not fix it in a certain

area for a long time. If you feel any fever or discomfort, please move

the irradiation slowly. (3) During the use of the device, collision and impact shall be prevented. (4) If you find that the brightness of the lamp is weak or not bright,

please contact our after-sales service. (5)Please, stay away from shortwave or microwave equipment while using the device. (6)Led irradiation may interfere with the work of the pacemaker. It is

Using the device during pregnancy is not recommended, unless otherwise advised by your physician. (7) Usage of the device by those with a pacemaker, implanted defibrillator, insulin pump is prohibited, because of possibly arising interference!

(8) The device shall not be disassembled by yourself, and large area of metal objects shall not be irradiated at close range. (9) Please do not use the device in the environment of explosive anesthetics or other in flammables.

(10) Please be sure to keep this product out of children's reach. (11) After use, a small number of people with thirst, dry eyes, rapid pulse, even dizziness and other symptoms should stop treatment.

Observe for a period of time, if the symptoms disappear, continue

treatment,if there is still,please stop.

exercises are necessary to enhance the metabolism.

Philippines, Singapore, Vietnam, Brazil, Colombia and Chile



sales guarantee services for all customers.



1)If you don't understand anything about our products, please I.INTRODUCTION 2.Benefits of the Cold Laser Therapy Device ———4

est return time limit of the shopping platform. Even if it
eeds the latest return time of the platform, we also provide
or-sales guarantee services, including return, refund.

1. Necessity of Cold Laser Therapy
2. Applicable Area
3. How Long Does it Take to Heal 3. How Long Does it Take to Heal -----




Solar therapy was first used by the Egyptians, and forms of l therapy were also practiced by the ancient Greeks, Chinese and

Red light resides on the visible end of the spectrum, which mean that you can see it. It is, however, on the "long end" of the whole spectrum in general.Red light waves are the longest wavelengths that are still visible to the human eye. The red light longer th

wavelength and the higher on the spectrum, the more effective light is at penetrating your body to a significant cellular level. Short red light penetrates your body to a depth of about 8 to 1 millimeters,reaching the mitochondria of your skin cells and othe

600 700 800 900

Wavelength in nm

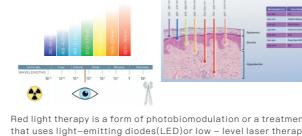

Is that are located closely beneath your skin to stimulat

Creating a temporary, low-dose metabolic stress (know works) that ultimately builds up the anti-inflammatory, antioxidant and cell defense systems of the cell.


1.Treatment Principle What's red light?

change and healing.

red light therapy and the most suitable wavelength is. Obviously,red light is best suited for treating skin conditions,

Then maybe someone will wonder that what we can benefit from Yet, the reason it can benefit all these radically different

It's important that almost all wavelengths of red LED light will energy. provide benefits,but wavelengths from 650 to 850nm demonstrate the most healing power.

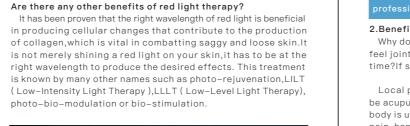


mbar spine, shoulders and neck, and soreness in the arms and

Secondly, the technology of red light therapy has been significantly knees. It is easy to cause diseases such as cervical and lumbar mature and has a very significant effect on local pain. spine. This is because of sitting for a long time for a long time, aintaining an improper posture, or because of work, it is easy cause frozen shoulder, mouse hand, lumbar disc herniation,

conditions is actually quite simple:The health of every organ and every cell in the body depends on the energy being produced b the mitochondria in those cells. Thus, because red/NIR light therapy

muscle strain, which will be accompanied by local soreness and pain.

2.Benefits of the Cold Laser Therapy Device Why do we need cold laser therapy in our daily life?Do you often feel joint pain, especially after sitting or working hard for a l time?If so,then local pain is damaging your body. be acupuncture-like pain. Pain itself is a signal. When the whole body is uncomfortable, it is conveyed through pain, such as knee

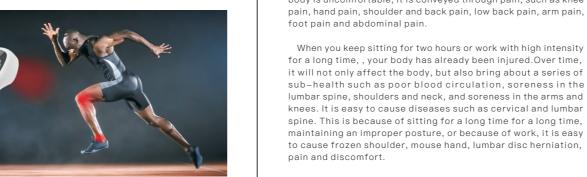

Question:Can you use vitamin C with red light therapy?

Answer: If your primary skincare goal is to brighten you

complexion with a boost of radiance, the brightening benefit

energy being produced by the mitochondria in those cells.those

cells will work better when the mitochondria are producing more

versatile addition to nearly any skincare routine.

Question:How did we get our start?

energy.

Question: Why do we love what we do?

Question: What makes our products unique?

of Vitamin C pair nicely with red light therapy. We could all use

a little extra glow from time to time –this makes Vitamin C a

Improve Joint Health and Combat Arthritis Improve Heart Health

Decrease PainIncrease Pancreas HealthImprove Immunity

which will produce pain under the stimulation of inflammation. 2.Excessive activity or high-intensity physical labor will cause

3. Frequent repetition of an action causes local muscle strain, which can cause to obvious pain. For example, frequent head lowering can lead to neck muscle spasm, which can lead to neck muscle soreness.

obvious symptoms of redness, swelling, heat and pain. There are many ways to relieve pain. First, rest, avoid excessive exertion. Reducing psychological tension through rest can alleviate

5. When bacterial infection occurs in local soft tissue, there will be

 Increase Fat Loss, Reduce Cellulite (And Burn Off Stubborn Fat) • Enhance Muscle Gain, Strength, Endurance, and Recovery Combat Inflammation (and Potentially Inflammation-Related Diseases)

Loco-motor diseases are the group of disorders which are t most common and affect the highest number of patients within the worldwide almost 30 % of the adult population suffers from locomotor pains, and one out of 5 persons is under long term

 No Adverse Short or Long Term Side Effects actual tasks of specialist care: consultation and special medical Easily Self-Administered in Your Own Home Long Term Application of the Device does not Have any Proven system.

as causes of death, even though this fact often remains hidden in mortality statistics.

supply of tissues is further degraded, the removal of metabolic

n the area of these types of treatments, and constitutes the perfec combination of the most effective therapies, lightweight, portable

1.Necessity of Cold Laser Therapy

Despite the fact that people's living standard and diet ha

a rather significant degree, they also represent a serious risk factor

and rechargeable device. Light irradiation using a low power density

any) side effects. Compared with specialized nursing in hospitals. portable handheld laser equipment is a very economical way of treatment.Cost performance is definitely more than you think,So you can easily relieve the pain without leaving home. Features of the Cold Laser Therapy Device as Below:

treatment, because of a rheumatologic disease and rheumatoid arthritis (joint inflammation). eumatology, to a lesser extent orthopedic and neurology specialist care providers, preventing the specialists from performing the

omfortably in pure physical therapy, and bring you a healthy and

dysfunctional cells differently than healthy cells. It even has th

applied in situations that can be apparently paradoxical because it

of laser therapy, which means that application to one site of the body

can produce an improvement of a condition in another distant body part that did not receive light. Systemic effects can be explained by ocal effects of light that can be transferred to other sites through


can sometimes be used to stimulate cells and tissues, and in ot

capacity to irradiate the blood (and affect things like inflammatory

Skin metabolism time: four to six months; Rone metabolism time: more than seven years

Products from one turn of the cycle Malate cycle α-Ketoglutarate Succinate (ATP)

From the above information, we can roughly understand why some diseases recover quickly, some more slowly. Therefore, in

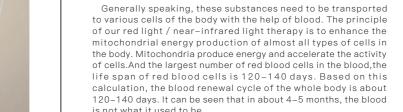

Similarly, the metabolic cycle of human tissues and organs is different. it affects the blood circulation channels of the whole body, and process needs to be maintained for a long time.


the curative effect.

3.Main Technical Parameters

1*User Manual 1*Portable Suitcase 1*USB Power Cable 1*Goggles 1*Charger Laser medium GaAlAs Semiconductor Laser wavelength 650nm+20nm/808nm+20nm Number of laser diodes 12pcs 808nm laser diodes Maximum output power 100, 150mW+20% r 650nm laser diode 5-10mW+20% no sotting = 5.00 minutes (5 minutes as as

NDHELD LIGHT THERAPY DEVICE

will be "ON", but it doesn't mean that it needs to be plugged n before it can be used. Our device has a Built-In battery inside, Type-C interface, supports USB charging, and can be used without plugging in.

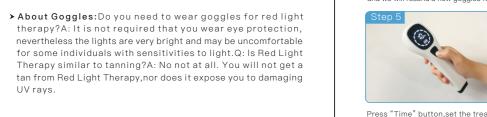



your problem within 24 hours. Our after-sales service is lifelong. It is not limited to the II.COLD LASER THERAPY DEVICE INTRODUCTION

Welcome to consult anytime!



III.MAIN STRUCTURAL COMPONENTS ------12

VI. TREATMENT RECOMMENDATION & GUIDELINES ----- 22

VIII.MAINTENANCE AND CLEANING -------27 IX.TRANSPORTATION AND STORAGE ------28

XI.CUSTOMER SUPPORT ------29

Question:Does red light therapy work for back pain? Answer: CONCLUSION. The IR therapy unit used was demonstrated to be effective in reducing chronic low back pain and no adverse

Question:ls red light therapy good for eyes?

Answer: Conclusion: Red Light Helps Protect and Heal Your Eyes

improvement to the look of your skin. It can be brighter clearer

and more even-toned. It can shrink your pores, reduce your

Answer: When the light is applied to your skin the epidermi absorbs it and then stimulates collagen proteins. In theory, mor

and blemishes, and tighten sagging skin, resulting in a more

Red light therapy treatments, with both red and near infra wavelengths, aren't just safe for eyes. They" re also a pro way to help heal eye injuries, reduce inflammation, and pro against vision loss, as shown in numerous peer-reviewed clin

Question:Does red light therapy work on face?

Question:What does light therapy do for skin?

Question:Does red light therapy tighten loose skin?

Question:Why choose our instrument?

Question: Why is this called cold laser? Does it have any heat

Sports Injuries & Sprain Skin Wound Healing Rheumatoid Arthritis

Reduce Joint Pain and Alleviate Arthrtis











Question: I want to use it on the highest power. How do you Question:How does red and near-infrared (NIR) light therapy

stimulating cell function. Even though cold lasers have a thermal effect and have an output of 1MW or above, but cold (or) Low-

Level Lasers do not have a thermal effect on Tissue, lasers that gentle heat but not cold.

Answer:Red/NIR light therapy work to enhance mitochondrial energy production in essentially every type of cell in the body, it can enhance the cellular processes and cellular health of potentially almost every type of cell in the body.

or is it completely cold?

Answer:Cold laser therapy is low-intensity laser therapy that

Answer:Infrared light therapy at home is just as effective to a sale or one which has a Red Light Thorapy going to a salon or spa which has a Red Light Therapy Device. wrinkles and give you a more youthful appearance.

> Question:Does it have any side-effects by using the therapy? Answer:In numerous studies,red and near infrared light has been

Question: What is the wavelength (nm) of this product? nswer:Perfect combination of 660nm & 850nm wavelengths The light therapy product uses the most effective wavelengt of visible red light (660nm) and Near infrared light (850nm) effects were observed. (Invisible) to help the pet skin and body recover.

Answer:-Pain Relief

milliwatts (ten – one thousandths of a watt), so you may feel a type of light.

with any side effects, at least if used short-term and as directed. This therapy is not toxic, not invasive and not as harsh as some topical skin treatments. Unlike the cancer-causing ultraviolet (UV) light from the sun or tanning booths, RLT doesn't use this described as the fountain of youth, red light therapy has been embrace new day. Light therapy device is the 100% natural

> Question: Is Red light therapy similar to tanning? Answer:No not at all. You will not get a tan from Red Light Therapy,

clinically shown to reduce wrinkles, promote the healing of scars solution for relieving the pain.

fer to the table below for detailed treatment guidelines:

TREATMENT GUIDELINES OF RED LIGHT DEVICE FOR PAIN RELIEF therapy device is the 100% natural solution for relieving the pain.

treatments per day with shorter duration. Twice/ 10days 2 3days ia 20 min Twice/per day 7days 2 2days treatments), in this case it is recommended to take a minimum

Close to the treatment site, change the angle of irradiation around 15 minutes long brake between treatments. 30 min Twice/per day 5days 4 2days 30 min Twice/per day 5days 4 2days r strain 30 min+ 1 wice/per day 3days 2 2days rain 30 min+ Twice/per day 5days 2 2days can be continued.At this period once a day and 30 minutes each time.Continuous treatment for 5 days. time.Continuous treatment for 5 days.

of the treatment also depends on the severity of the pain. time.Continuous treatment for 5 days. the problem ceases. Long term application does not have any

Physical therapy is a long-term effective painless and sideproven side-effects.For normal wound ,use this device twice effect-free treatment, especially for people with limited mobility, Treatment Site: so that you can easily relieve pain at home. Nothing beats a Close to the treatment site, change the angle of irradiation strong constitution for being ready to embrace new day.Light around the lesion every 10 minutes of treatment.

s 3 0 min+ Twice/per day 7 days 3 2 days 10 minutes and a maximum of 30 minutes. The number of daily each time. Continuous treatment for 7 days. treatments is a minimum of 3. In more serious cases several treatments can be performed in succession (as many as 10 Treatment site:

> (1) Treatments for Acute Soft Tissue Injury
> Treatment instruction: st twice a day and about 30 minutes each time.10 days is one treatment session and3 days' interval between 2 sessions. It is recommended to use 2 treatment sessions. After the curative effect appears, the consolidation treatment be continued. At this period once a day and 30 minutes each

(2)Treatments for Shoulder Periarthritis and Knee Arthritis We suggest twice a day and about 30 minutes each time. 7 days is one treatment session and 2 days' interval between 2 sessions. It is recommended to use 3 treatment sessions.

After the curative effect appears, the consolidation treatment

each time.Continuous treatment for 3 days. Treatment site:close to the treatment site.

various problems and their severity. At the same time, the effect can be continued. At this period once a day and 30 minutes each insist on using it continuously when you want it to have a

the lesion every 10 minutes of treatment.

After the curative effect appears, the consolidation treatmen

The duration of treatments may be different, depending on the

After the curative effect appears, the consolidation treatment

• The therapy effect varies from person, to person, Please significant therapeutic effect on pain relief. The short-term effect may not be obvious that depends on your body condition. If you have any questions please feel free to contact us.

the curative effect appears, the consolidation treatment can

(4)Treatments for Cervical and Lumbar Osteoarthritis eatment instruction: le suggest twice a day and about 30 minutes each time. 5

Treatment site: Close to the treatment site.

 After switching it off always wipe the device with a moistened. cotton wiping cloth,thus preserving hygiene.Do not use

and do not expose it to direct sunlight, because that reduce the device's lifespan. ●Please,do not place hard or sharp objects on or near the

emitting ports or LED light probes.Do not immerse the device

(1)The host can be replaced for free with warranty card in the case of quality problems under normal use within one year device, which may damage its surface or the laser treatment

> s immediately, and we can reissue it free for you (2)Damages caused by improper operation or mistakes. (3)Breakdowns caused by self-disassembling. The device's repair was not performed by the professionals of the service or the device was dismantled outside of the service.

your understanding and support!

other accessories.

Email Format Template:

Laser therapy involves red light, infrared light and low–intensity laser into the human body by irradiation. Product Front View Goggles USB Power Cable Wipes Cloth Charger Therefore, when the red light therapy directly acts on the cells

Product Front View Product Back View

(2)The company undertakes the device's repair after the expiration of the warranty period as well, Including but not

activities scratches, soaking or unauthorized replacement of

(5)Damages caused by natural disasters.(6)The damage of shell,manual,packaging,warranty card and

During use you cannot sense

The skin is too dry

Moisten the skin with a wet cloth, or use gel on the affected skin surface Stimulation is weak is not properly in contact with the skin head, so it will contact the affected skin surface with the largest surface

Out of power Contact the distributor for Other malfunction Guarantee, repair or replacement

Lamp Flicker Environmental interference

The light suddenly went off

Low input voltage Check for flooding

Check that the power supply is plugged in properly.

Cut of power Recharge pow

of the light Low input voltage Check compatibility,replace Incorrect contact

LCD display turns
off during use

Normal program

No operation, the LCD display will go out after 2 minutes, it's normal.

No operation, the LCD display will go out after 2 minutes, it's normal.